


# AD-Line™



Advanced tubing line with multiple options.

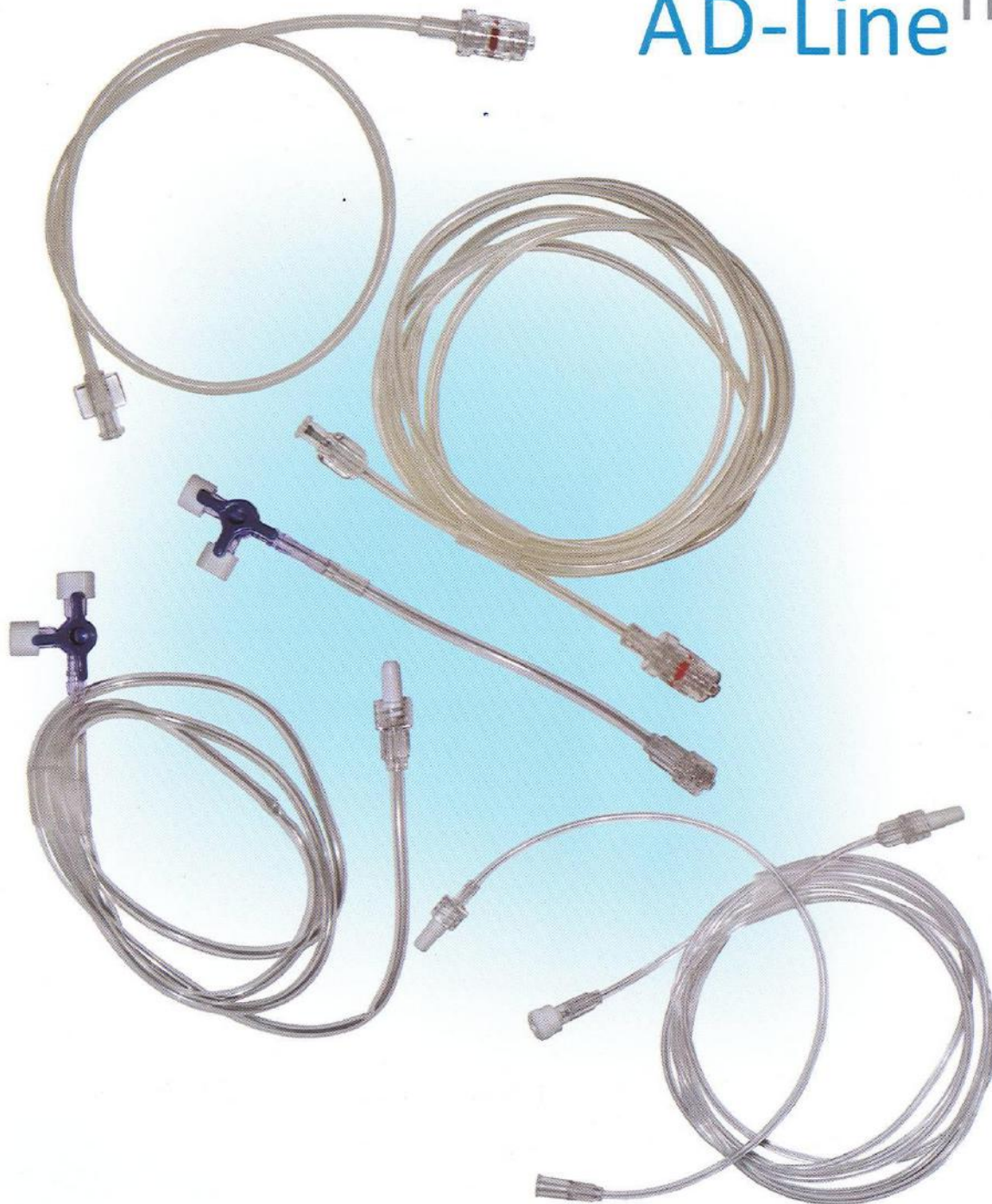


Changing lives for the better...

## AD-LINE™ HP 1200 PSI CONNECTING TUBING 1200 PSI High Pressure Injection Tubing

### Description:

- Braided tubing
- Used for power injection applicator
- Combines high pressure capacity with superior flexibility
- Available with or without high pressure rotating adapter
- Latex free
- The following High pressure line are available in Male / Female Luer Lock.

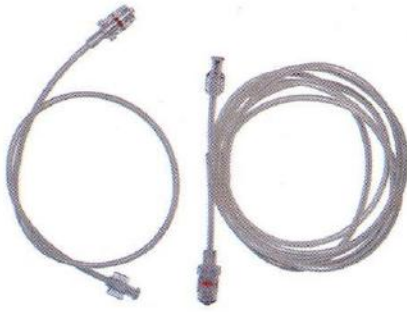

| Catalog No | Description          |
|------------|----------------------|
| ADLH25MF   | 25cm Male/Female LL  |
| ADLH50MF   | 50cm Male/Female LL  |
| ADLH100MF  | 100cm Male/Female LL |
| ADLH120MF  | 120cm Male/Female LL |
| ADLH150MF  | 150cm Male/Female LL |
| ADLH200MF  | 200cm Male/Female LL |

## AD-LINE™ 3 Way Stop Cock with Extension Line

### Features:

- Superior quality 3-Way Stop Cock now available with integrated extension tubing.
- Extension tube provides separation of administration and insertion sites.
  - i) Minimizes mechanical irritation of blood vessel being used for access.
  - ii) Reduces thrombotic events such as thrombophlebitis.
  - iii) Increases patient comfort.
- 3-Way Stop Cock facilitates multiple lines through single IV access.
- Ideal for administering contrast medium/fluids during imaging procedures that require patient to enter a gantry; e.g. CT and MRI Scans.
- High quality kink-resistant tubing with smooth inner surfaces for laminar, non-turbulent flow.
- Available in various tube lengths based on customer requirements



### 3 Way Stop Cock with Extension line

| Catalog No | Length cm |
|------------|-----------|
| AL3W10MF   | 10 cm     |
| AL3W25MF   | 25 cm     |
| AL3W50MF   | 50 cm     |
| AL3W100MF  | 100 cm    |
| AL3W120MF  | 120 cm    |
| AL3W150MF  | 150 cm    |
| AL3W200MF  | 200 cm    |


## AD-LINE™ Pressure Monitoring Line

Pressure monitoring tubing is non-expanding PVC tubing with a rigid wall thickness which will ensure high quality wave forms with improved clarity. Available with both male-to-male & male-to-female connectors from 15 cm to 200 cm in length.



| Catalog No | Length cm | Tubing ID | Tubing OD | Proximal Hub Style |
|------------|-----------|-----------|-----------|--------------------|
| ADLP15MF   | 15 cm     | 1.0 mm    | 2.8 mm    | Male/Female        |
| ADLP25MF   | 25 cm     | 1.0 mm    | 2.8 mm    | Male/Female        |
| ADLP50MF   | 50 cm     | 1.0 mm    | 2.8 mm    | Male/Female        |
| ADLP100MF  | 100 cm    | 1.0 mm    | 2.8 mm    | Male/Female        |
| ADLP120MF  | 120 cm    | 1.0 mm    | 2.8 mm    | Male/Female        |
| ADLP150MF  | 150 cm    | 1.0 mm    | 2.8 mm    | Male/Female        |
| ADLP200MF  | 200 cm    | 1.0 mm    | 2.8 mm    | Male/Female        |

Not for sale in USA

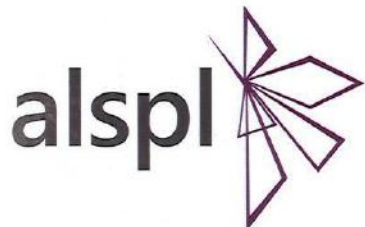

Changing lives for the better...

### ADVANCED LIFE SCIENCES PVT. LTD.

D-22, Okhla Industrial Area, Phase-I, New Delhi-110020, INDIA

Tel. : +91 (0)11 - 26812108, 26371127, 26812107

Fax : +91 (0)11 - 26371128

Email: customerservice@alspl.com; Website: [www.alspl.com](http://www.alspl.com)