

ARROW



Arrow

Super Arrow-Flex Sheath Introducer

Flexibility you can see. Strength you can feel.



Teleflex

Complete vascular access without kinking or collapsing

The Flexible, Kink-Resistant Super Arrow-Flex Sheath Introducer

Arrow Sheath Introducers from Teleflex continue to be the choice of many clinicians who require flexibility in sheath design in many areas of the world. The Super Arrow-Flex Sheath Introducer gives you complete vascular access, thanks to its wire-reinforced design and construction. From easy insertion and handling to the highly visible radiopaque tip*, Super Arrow-Flex Sheath Introducers are designed to deliver optimal performance from start to finish. And with a wide range of sheath sizes and lengths, Super Arrow-Flex Sheath Introducers give you flexibility of choice for multiple vascular applications.

Effectiveness even in the most challenging cases

- Exclusive coil-wire design allows the sheath to flex at any point in any direction without kinking or losing support
- Excellent steerability helps negotiate tortuous anatomies
- Tapered dilator provides pushability to the contralateral leg and Superficial Femoral Artery (SFA)

Easy placement and insertion

- Radiopaque tip marker enhances visualization during insertion and enables precise deployment of a catheter just beyond the sheath tip.*
- Highly radiopaque sheath, yet spacing between coils enables visualization through the sheath

- Smooth transition between the sheath and dilator helps to prevent peel-back
- Hydrophilic surface helps facilitate smooth insertion
- Valve lubrication eases catheter placement through the sheath


Range of products for interventional procedures

- Specific designs suitable for renal, carotid and transseptal procedures
- Wide range of products providing reliable access for cross-over, TIPS and antegrade brachial procedures
- Variety of product options from 5 Fr. To 10 Fr. and 11 cm to 90 cm
- Sheath available with and without guide wire
- Color-coded for convenience

Safety and peace of mind

- Ability to be used as a guiding sheath may eliminate the need for a guiding catheter (includes obturator cap)
- Color-coded hubs for easy identification of sheath sizes
- Tri-cusped hemostasis valve facilitates catheter manipulation with minimal bleed-back
- Sheath and dilator together to enable simultaneous, secure insertion

*In all 35 cm and longer except 5 Fr.

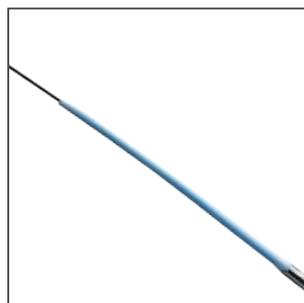

Each CP set includes:

- One Super Arrow-Flex Sheath with Integral Side Port/Hemostasis Valve and attached 3-Way Stopcock
- One 17 3/4" (45 cm) long dual purpose "J" tip spring-wire guide in holder with tip straightener. Short-length sheath set (P) has .021 dia. and all others have .035" dia. wire guide
- One Vessel Dilator with SnapLock Feature
- One Obturator Cap

Each CL set includes::

- One Super Arrow-Flex Sheath with Integral Side Port/Hemostasis Valve and attached 3-Way Stopcock
- One Vessel Dilator with SnapLock Feature
- One Obturator Cap

Super Arrow-Flex Sheath Introducer Sets with Wire Guides

COLOR CODE	PRODUCT NO.	SHEATH SIZE	SHEATH LENGTH	SETS PER CASE
●	CP-07511	5 Fr.	11 cm	10
●	CP-07611	6 Fr.	11 cm	10
●	CP-07711	7 Fr.	11 cm	10

Super Arrow-Flex Sheath Introducer Sets Without Wire Guides

COLOR CODE	PRODUCT NO.	SHEATH SIZE	SHEATH LENGTH	SETS PER CASE
6 Fr. Introducers				
●	CL-07624	6 Fr.	24 cm	5
●	CL-07645	6 Fr.	45 cm	5
●	CL-07665	6 Fr.	65 cm	5
●	CL-07690	6 Fr.	90 cm	2
7 Fr. Introducers				
●	CL-07724	7 Fr.	24 cm	5
●	CL-07745	7 Fr.	45 cm	5
●	CL-07765	7 Fr.	65 cm	5
●	CL-07780	7 Fr.	80 cm	2
●	CL-07790	7 Fr.	90 cm (Carotid)	2
8 Fr. Introducers				
●	CL-07824	8 Fr.	24 cm	5
●	CL-07845	8 Fr.	45 cm	5
●	CL-07865	8 Fr.	65 cm	5
●	CL-07880	8 Fr.	80 cm	2
9 Fr. Introducers				
●	CL-07980	9 Fr.	80 cm	2
10 Fr. Introducers				
●	CL-07080	10 Fr.	80 cm	2



Teleflex is a global provider of medical technologies designed to improve the health and quality of people's lives. We apply purpose driven innovation – a relentless pursuit of identifying unmet clinical needs – to benefit patients and healthcare providers. Our portfolio is diverse, with solutions in the fields of vascular and interventional access, surgical, anesthesia, cardiac care, urology, emergency medicine and respiratory care. Teleflex employees worldwide are united in the understanding that what we do every day makes a difference. For more information, please visit teleflex.com.

Teleflex is the home of Arrow, Deknatel, Hudson RCI, LMA, Pilling, Rüsçh and Weck – trusted brands united by a common sense of purpose.

Corporate Office

Phone +1 610 225 6800, 550 E. Swedesford Road, Suite 400, Wayne, PA 19087, USA

Regional Offices

United States: Phone +1 919 544 8000, Toll Free 866 246 6990, cs@teleflex.com, 3015 Carrington Mill Boulevard, Morrisville, NC 27560, USA

Latin America: Phone +1 919 433 4999, la.cs@teleflex.com, 3015 Carrington Mill Boulevard, Morrisville, NC 27560, USA

International: Phone +353 (0)9 06 46 08 00, orders.intl@teleflex.com, Teleflex Medical Europe Ltd., IDA Business and Technology Park, Dublin Road, Athlone, Co Westmeath, Ireland

Australia 1300 360 226

Austria +43 (0)1 402 47 72

Belgium +32 (0)2 333 24 60

Canada +1 (0)800 387 9699

China (Shanghai) +86 (0)21 6163 0965

China (Beijing) +86 (0)10 6418 5699

Czech Republic +420 (0)495 759 111

France +33 (0)5 62 18 79 40

Germany +49 (0)7151 406 0

Greece +30 210 67 77 717

India +91 (0)44 2836 5040

Italy +39 0362 58911

Japan +81 (0)3 6632 3600

Korea +82 2 536 7550

Mexico +52 55 5002 3500

Netherlands +31 (0)88 00 215 00

New Zealand 0800 601 100

Portugal +351 22 541 90 85

Singapore (SEA non-direct sales countries) +65 6439 3000

Slovak Republic +421 (0)3377 254 28

South Africa +27 (0)11 807 4887

Spain +34 918 300 451

Switzerland +41 (0)31 818 40 90

United Kingdom +44 (0)1494 53 27 61

For more information, please visit teleflex.com.

Teleflex, the Teleflex logo, Arrow, Deknatel, Hudson RCI, LMA, Pilling, Rüsçh, SnapLock, Super Arrow-Flex and Weck are trademarks or registered trademarks of Teleflex Incorporated or its affiliates, in the U.S. and/or other countries. Information in this document is not a substitute for the product Instructions for Use. The products in this catalogue may not be available in all countries. Please contact your local representative. All data current at time of printing (05/2017). Subject to technical changes without further notice. © 2017 Teleflex Incorporated. All rights reserved.