



SUPER ARROW-FLEX[®] SHEATH INTRODUCER

Flexibility you can see. Strength you can feel.


COMPLETE VASCULAR ACCESS WITHOUT KINKING OR COLLAPSING

THE FLEXIBLE, KINK-RESISTANT SUPER ARROW-FLEX SHEATH INTRODUCER

The ARROW brand of Teleflex sheath introducers continues to be the choice of many clinicians who require flexibility in sheath design in many areas of the world. Super Arrow-Flex gives you complete vascular access, thanks to its wire-reinforced design and construction. From easy insertion and handling to the highly visible radiopaque tip*, Super Arrow-Flex is designed to deliver optimal performance from start to finish. And with a wide range of sheath sizes and lengths, Super Arrow-Flex gives you flexibility of choice for multiple vascular applications.

RANGE OF PRODUCTS FOR INTERVENTIONAL PROCEDURES

- specific designs suitable for renal, carotid, and transseptal procedures
- wide range of products provides reliable access during crossover, TIPS, and antegrade brachial procedures
- variety of product options from 5Fr. to 11Fr. and 7.5 cm to 100 cm
- sheath available with and without guide wire
- colour coded for convenience

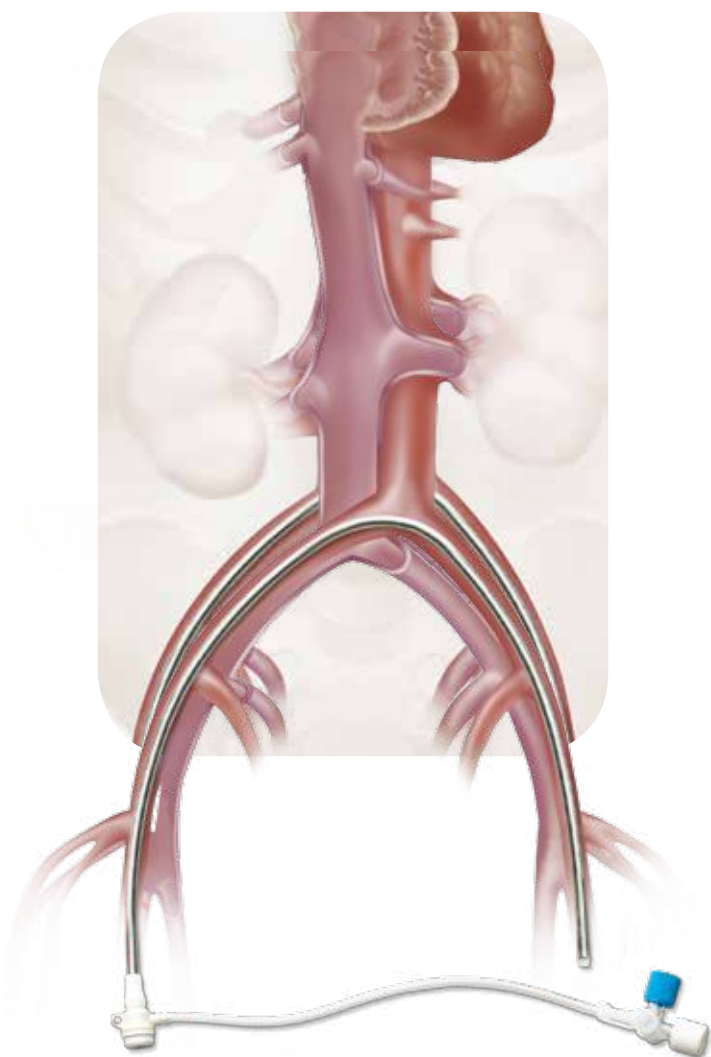

*35 cm and longer, except 5 Fr.

MAXIMUM EFFECTIVENESS EVEN IN THE MOST CHALLENGING CASES

- exclusive coil-wire design allows the sheath to flex at any point and in any direction without kinking or losing support
- excellent steerability helps negotiate tortuous anatomies
- tapered dilator provides exceptional pushability to the contralateral leg and the superficial femoral artery (SFA)
- increased patient comfort and reduced back pain post-procedures due to the flexible design¹



¹ Waksman et al. Randomized comparison of flexible versus nonflexible femoral sheaths on patient comfort after angioplasty. *American Heart Journal*. 131:6, 1076-1078.

SUPER ARROW-FLEX SHEATH INTRODUCER SETS WITH WIRE GUIDES

5-10 FR. INTRODUCERS

REF.	SHEATH SIZE	SHEATH LENGTH	COLOUR CODED HUB	SETS/CASE
CP-07511-P	5 Fr.	7.5 cm	grey	10
CP-07511	5 Fr.	11 cm	grey	10
CP-07611-P	6 Fr.	7.5 cm	green	10
CP-07611	6 Fr.	11 cm	green	10
CP-07711	7 Fr.	11 cm	orange	10
CP-07811	8 Fr.	11 cm	blue	10
CP-07911	9 Fr.	11 cm	black	10
CP-07011	10 Fr.	11 cm	white	10

EACH "CP" SET INCLUDES:

- Super Arrow-Flex sheath with integral side port/haemostasis valve and attached 3-way stopcock
- 17 3/4" (45 cm) long dual purpose "J" tip spring-wire guide in holder with tip straightener. Short-length sheath set (P) has .021"-dia. and all others have .035"-dia. wire guide
- vessel dilator with snap-lock feature
- obturator cap

SUPER ARROW-FLEX SHEATH INTRODUCER SETS WITHOUT WIRE GUIDES

5 FR. INTRODUCERS

REF.	SHEATH SIZE	SHEATH LENGTH	COLOUR CODED HUB	SETS/CASE
CL-07511	5 Fr.	11 cm		10
CL-07524	5 Fr.	24 cm		5
CL-07545	5 Fr.	45 cm	grey	5
CL-07565	5 Fr.	65 cm		5
CL-07590	5 Fr.	90 cm		2

6 FR. INTRODUCERS

REF.	SHEATH SIZE	SHEATH LENGTH	COLOUR CODED HUB	SETS/CASE
CL-07611	6 Fr.	11 cm		10
CL-07624	6 Fr.	24 cm		5
CL-07635	6 Fr.	35 cm		5
CL-07645	6 Fr.	45 cm	green	5
CR-07645*	6 Fr.	45 cm (renal)		5
CL-07665	6 Fr.	65 cm		5
CL-07690	6 Fr.	90 cm		2

7 FR. INTRODUCERS

REF.	SHEATH SIZE	SHEATH LENGTH	COLOUR CODED HUB	SETS/CASE
CL-07711	7 Fr.	11 cm		10
CL-07724	7 Fr.	24 cm		5
CL-07735	7 Fr.	35 cm		5
CL-07745	7 Fr.	45 cm		5
CR-07745*	7 Fr.	45 cm (renal)		5
CR-07745-NT**	7 Fr.	45 cm (renal)	orange	5
CL-07765	7 Fr.	65 cm		5
CL-07780	7 Fr.	80 cm		2
CL-07790-R*	7 Fr.	90 cm (carotid)		2
CL-07700	7 Fr.	100 cm		2



EASY PLACEMENT AND INSERTION

- radiopaque tip marker enhances visualisation during insertion and enables precise deployment of a catheter just beyond the sheath tip
- highly radiopaque sheath, yet spacing between coils enables visualisation
- smooth transition between sheath and dilator helps prevent peel-back
- hydrophilic surface facilitates smooth insertion
- valve lubrication eases catheter placement through the sheath



SAFETY AND PEACE OF MIND

- capability to be used as a guiding sheath may eliminate the need for a guiding catheter (includes obturator cap)
- colour-coded hubs for easy identification of sheath sizes
- tri-cuspid haemostasis valve facilitates catheter manipulation with minimal bleed-back
- sheath and dilator lock together to enable simultaneous and secure insertion



8 FR. INTRODUCERS

REF.	SHEATH SIZE	SHEATH LENGTH	COLOUR CODED HUB	SETS/ CASE
CL-07811	8 Fr.	11 cm		10
CL-07824	8 Fr.	24 cm		5
CL-07835	8 Fr.	35 cm		5
CL-07845	8 Fr.	45 cm		5
CT-07860	8 Fr.	60 cm (transseptal)	blue	2
CL-07865	8 Fr.	65 cm		5
CL-07880	8 Fr.	80 cm		2
CL-07890-R*/***	8 Fr.	90 cm (carotid)		2
CL-07800	8 Fr.	100 cm		2

9 FR. INTRODUCERS

REF.	SHEATH SIZE	SHEATH LENGTH	COLOUR CODED HUB	SETS/ CASE
CL-07911	9 Fr.	11 cm		10
CL-07924	9 Fr.	24 cm		5
CL-07965	9 Fr.	65 cm	black	5
CL-07980	9 Fr.	80 cm		2
CL-07900	9 Fr.	100 cm		2

10 FR. INTRODUCERS

REF.	SHEATH SIZE	SHEATH LENGTH	COLOUR CODED HUB	SETS/ CASE
CL-07011	10 Fr.	11 cm		10
CL-07024	10 Fr.	24 cm		5
CL-07035	10 Fr.	35 cm	white	5
CL-07045	10 Fr.	45 cm		5
CL-07065	10 Fr.	65 cm		5
CL-07080	10 Fr.	80 cm		2

11 FR. INTRODUCERS

REF.	SHEATH SIZE	SHEATH LENGTH	COLOUR CODED HUB	SETS/ CASE
CL-71165	11 Fr.	65 cm		5
CL-71180	11 Fr.	80 cm	yellow	2

EACH "CL" SET INCLUDES:

- Super Arrow-Flex sheath with integral side port/ haemostasis valve and attached 3-way stopcock
- vessel dilator with snap-lock feature
- obturator cap

* These products include a special Tuohy-Borst "Y" adapter with side arm extension.

** Tuohy-Borst "Y" adapter is not included with this product.

*** Contact your ARROW representative to check availability.

Teleflex is a leading global provider of specialty medical devices used for diagnostic and therapeutic procedures in critical care, urology and surgery. Our mission is to provide solutions that enable healthcare providers to improve outcomes and enhance patient and provider safety. We specialise in devices for general and regional anaesthesia, cardiac care, respiratory care, urology, vascular access and surgery and we serve healthcare providers in more than 130 countries. Teleflex also provides specialty products for medical device manufacturers.

Our well known brands include ARROW[®], DEKNATEL[®], GIBECK[®], HUDSON RCI[®], KMEDIC[®], LMA[™], PILLING[®], PLEUR-EVAC[®], RÜSCH[®], SHERIDAN[®], TAUT[®], TFX OEM[®], VASONOVA[™] and WECK[®], all of which are trademarks or registered trademarks of Teleflex Incorporated.

Teleflex global operations: Austria, Belgium, Canada, China, Czech Republic, France, Germany, Greece, India, Ireland, Italy, Japan, Malaysia, Mexico, Netherlands, Portugal, Singapore, Slovak Republic, South Africa, Spain, Switzerland, United Kingdom, Uruguay and USA.

YOUR INTERNATIONAL CONTACTS:

TELEFLEX HEADQUARTER INTERNATIONAL, IRELAND

Teleflex Medical Europe Ltd., IDA Business and Technology Park,
Dublin Road, Athlone, Co Westmeath, Ireland
Phone +353 (0)9 06 46 08 00 · Fax +353 (0)14 37 07 73
orders.intl@teleflex.com

AUSTRIA +43 (0)1 402 47 72
BELGIUM +32 (0)2 333 24 60
CHINA (SHANGHAI) +86 (0)21 6381 6399
CHINA (BEIJING) +86 (0)10 6418 5699
CZECH REPUBLIC +420 (0)495 759 111
FRANCE +33 (0)5 62 18 79 40
GERMANY +49 (0)7151 406 0
GREECE +30 210 67 77 717
INDIA +91 (0)44-2836 5040
ITALY +39 0362 58 911
JAPAN +81 (0)3 3379 1511
NETHERLANDS +31 (0)88 00 215 00
PORTUGAL +351 22 541 90 85
SINGAPORE + 65 6447 5022
SLOVAK REPUBLIC +421 (0)3377 254 28
SOUTH AFRICA +27 (0)11 807 4887
SPAIN +34 918 300 451
SWITZERLAND +41 (0)31 818 40 90
UNITED KINGDOM +44 (0)1494 53 27 61

For detailed information see www.teleflex.com

The products in this catalogue may not be available in all countries.
Please contact your local representative. All data current at time of printing (02/2013).
Subject to technical changes without further notice.

94 10 43 - 00 00 01 · REV A · MC / WM · 02 13 02